

AEROTRAK® PORTABLE PARTICLE COUNTER MODEL 9500



The TSI AeroTrak® 9500 Portable Particle Counter offers best-in-class features and versatility. Perform ISO 14644-1 classifications and complete certifications of EU GMP facilities by sampling 1 m³ of air in just 10 minutes! Particle count data can be conveniently viewed on screen, downloaded using TrakPro™ Lite Secure software, or printed directly to the integrated printer. Compatible with TSI's accurate ventilation probes that measure air velocity/temperature/RH, the Model 9500 is a complete solution for your certification needs. The combination of an optional electronic filter scanner with up to 26 ft of tubing, makes this the best instrument for in-place filter testing.

This particle counter can generate pass/fail reports for ISO 14644-1, EU GMP Annex 1, and FS209E, via the onboard printer and the supplied TrakPro Lite Secure download software. Multiple particle counter configurations can be conveniently stored and uploaded as needed. It can be used as a stand-alone particle counter or integrated into a facility monitoring system like TSI's FMS 5.


The Model 9500 complies with all the stringent requirements set forth in ISO 21501-4. It is calibrated with NIST traceable PSL spheres using TSI's world-class Classifier and Condensation Particle Counters, the recognized standard for particle measurements. Backed with a standard two-year warranty and TSI's long-standing reputation for high quality, there are no other particle counters like it on the market today.

Features and Benefits

- + 0.5 to 25 µm size range
- + Up to six channels of simultaneous data
- + 3.53 CFM (100 L/min) flow rate
- + Complies with all requirements of ISO 21501-4
- + Best choice for filter testing in cleanrooms
- + Data storage: 250 zones, 999 locations, and 10,000 samples
- + USB, USB Storage Device and Ethernet data transfer
- + Industry's loudest audible alarm
- + Easy icon-driven touch screen operation
- + Onboard reports for ISO-14644-1, EU GMP Annex 1, and FS209E
- + External certification reports via TrakPro Lite Secure
- + Allows storage and recall of specific Recipes and Zone configurations
- + Displays up to three environmental parameters
- + Stainless steel enclosure
- + Lightweight
- + Stand-alone or integrate into a facility monitoring system



UNDERSTANDING, ACCELERATED

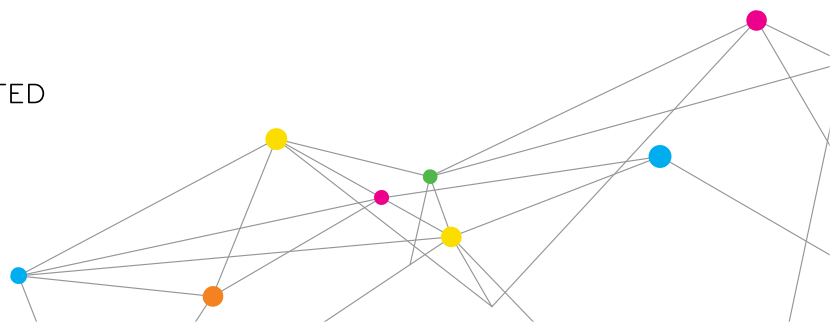

SPECIFICATIONS

AEROTRAK PORTABLE PARTICLE COUNTER MODEL 9500

9500		
Size Range	0.5 to 25 µm	250 Zones
Particle Channel Sizes	0.5, 0.7, 1.0, 3.0, 5.0, 10.0 µm (additional channel sizes available)	999 Locations
Size Resolution	<15% @ 0.5 µm (per ISO 21501-4)	10,000 sample records including: date, time, six total particle size channels, flow status, and instrument status
Counting Efficiency	50% at 0.5 µm, 100% for particles > 0.75 µm (per ISO 21501-4 and JIS)	Data Transfer
Concentration Limit	425,000 particles/ft ³ (15,000,000/m ³) @ 10% coincidence loss	Via USB storage device, connection to the TrakPro Lite Secure software through Ethernet or USB cable, or optional TSI FMS Software
Light Source	Long life laser diode	Reports
Zero Count	<1 count per 5 minutes (per ISO 21501-4 and JIS B9921)	Provides Pass/Fail on ISO 14644-1, EU GMP, and FS209E reports
Flow Rate	3.53 CFM (100 L/min) ±5% accuracy (meets ISO 21501-4 and JIS B9921)	Calibration
Flow Rate Control	Electronic, automatic closed loop (patented* flow control technology)	NIST traceable using TSI calibration system
Sample Tube Extension	Up to 26.2 ft (8 m)	Calibration Frequency
Sampling Modes	Manual, automatic, beep; cumulative/differential; count or concentration	Recommended minimum of once per year
Sampling Time	1 second to 99 hours	Dimension (H x W x D)
Sampling Frequency	1 to 9999 cycles or continuous	9.5 in. x 9.0 in. x 9.2 in. (24.1 cm x 23.0 cm x 23.2 cm) (without isokinetic inlet)
Exhaust	Internally HEPA filter	Weight
Vacuum Source	Internal pump	14.2 lbs (6.44 kg) with two batteries
Communication Mode	Modbus® TCP over Ethernet or USB	Power
Audible Alarm	Built-in; >85 dB at 1 meter (adjustable)	110 to 240 VAC universal power supply
Status Indicators	Low battery, flow, laser	Battery
Alarm Output	Dry contacts, closed when alarm is engaged	Removable/rechargeable Li-Ion batteries (2 battery slots)
Alarm Limits	Programmable for all particle channels	Battery Life
Display	VGA 5.7-inch (14.5 cm) touch-screen display	>2 hours of continuous operation
Languages	English, German, French, Spanish, Japanese, Chinese (simplified), Italian	Battery Recharge
Printer	Built-in thermal printer (also available without printer)	4 hours (when idle)
Printer Output	Prints in all available languages	External Surface
Environmental Sensor Interface	Supports TSI Probe Models: 960, 962, 964, 966 (air velocity/temp/RH) and Temp/RH probe	Stainless steel
External Chemical Resistance	Isopropyl alcohol, chlorinated solution, hydrogen peroxide	Operating Range
Software	Compatible with TrakPro Lite Secure and FMS 5 software	36° to 95°F (2° to 35°C), 20% to 95% RH noncondensing
Unit ID	Configurable IP address	Storage Range
Security	2-level password protection to lock out usage and configuration	32° to 122°F (0° to 50°C), up to 98% RH noncondensing
		Standards
		ISO 21501-4, CE, JIS B9921
		Warranty
		Two years, extended warranties available
		Included Accessories
		Printed QuickStart guide, operating manual on CD, power supply, battery, barbed inlet, isokinetic probe, tubing, purge filter, printer paper, USB cable and TrakPro Lite Secure Software
		Optional Accessories
		Electronic filter scanning probe, basic filter scanning probe, TSI velocity probes, Temp/RH probe, additional battery, dual port external battery charger, isokinetic probes, sample tubing, carrying case

*Patent Number 6,167,107

Specifications are subject to change without notice.

AeroTrak, TSI and the TSI logo are registered trademarks, and TrakPro is a trademark of TSI Incorporated.

Modbus is a registered trademark of Modicon, Inc.



UNDERSTANDING, ACCELERATED

TSI Incorporated - Visit our website www.tsi.com for more information.

USA Tel: +1 800 874 2811 **India** Tel: +91 80 67877200
UK Tel: +44 149 4 459200 **China** Tel: +86 10 8219 7688
France Tel: +33 1 41 19 21 99 **Singapore** Tel: +65 6595 6388
Germany Tel: +49 241 523030